

Patient Instructions for ER900 (AT)

To Record: Events to be recorded will be described by the physician.



1. Press the RECORD/SEND button for **two seconds** until an audible tone is heard, then release.

Notes:

- A. There are two RECORD/SEND buttons, you only need to push one of them for two seconds to activate the Monitor.
 - B. To stop a recording, press the RECORD/SEND button again for **two seconds**.
 - C. "RECORDING" will flash in the upper left hand corner of the display.
2. Hold as still as possible during recording but continue breathing.
 3. The recording is complete and ready to send when a phone ring is heard from the Monitor.

Notes:

- A. The display will also show that an event is stored in the Monitor. If the memory is full, follow instructions To Send and Erase Events

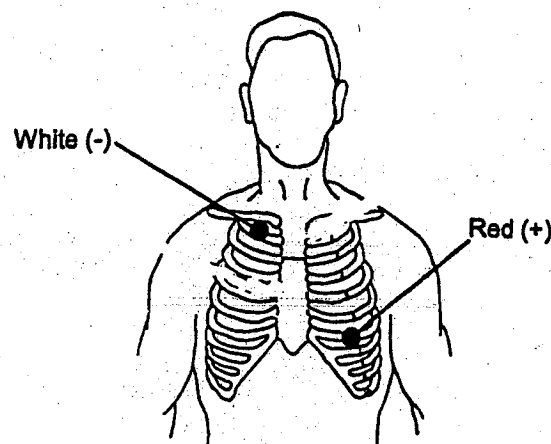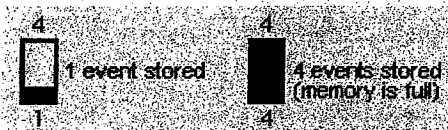

Looping Event Recorders
Common Electrode Placement

To Send:

1. Remove the Patient Cable from the Monitor.
2. Set the Monitor on a flat surface with the speaker hole up.
3. Call the receiving center.
 - A. Cell phones and VOIP phones **do not** work for the transmission.
4. Follow receiving center instructions.
 - A. When instructed, place the telephone **mouthpiece** over the Monitor speaker hole.



- B. Push the RECORD/SEND button for two seconds. The Monitor will make audible tones and "SENDING" will be shown in the upper left hand corner of the display while the recording is being sent.
5. When instructed that the events have been sent successfully, it is OK to erase events by reinserting Patient Cable. (The Monitor will restart and you will then hear a rising tone.)

WHEN YOUR MONITORING TIME IS COMPLETED, RETURN YOUR UNIT IN THE POSTAGE PAID ENVELOPE DIRECTLY TO LSSI. DO NOT RETURN THE UNIT TO THE CLINIC UNLESS OTHERWISE SPECIFIED BY YOUR PHYSICIAN OR CLINIC.

*****QUESTIONS/PROBLEMS PLEASE CONTACT LSSI @ 1-800-659-8151. TECHNICIANS AND NURSES ARE AVAILABLE 24/7/365 TO ASSIST YOU IN ALL YOUR NEEDS*****



(800) 659-8151

www.lifesupportinc.com